

Reference these catalogue numbers to order:


EMVL-100: Environmental Monitoring Vision Lids, 100 per case.

Part number	Name	Shelf Life	Storage Conditions
SKU-EMVL-100	EM Vision Lids	2 Years	Room Temperature

Product Description: The environmental monitoring vision lid is a transparent lid placed on top of the ET80 or ET80HT cassette after sampling is complete. It is used to protect the sampling surface during loading and testing on the Growth Direct® System. The vision lid should be firmly engaged to the media cassette. The vision lid allows images to be taken by the Growth Direct® System of the membrane on the enclosed media surface.

Manufacture: The vision lid product is manufactured in an ISO certified facility. The sterilization treatment is performed in an ISO certified facility.

Material: Vision Lid is made from Zeonor 1060R and the color is clear. The vision lids are packaged in 10 stacks of 10, in three Tyvek bags.

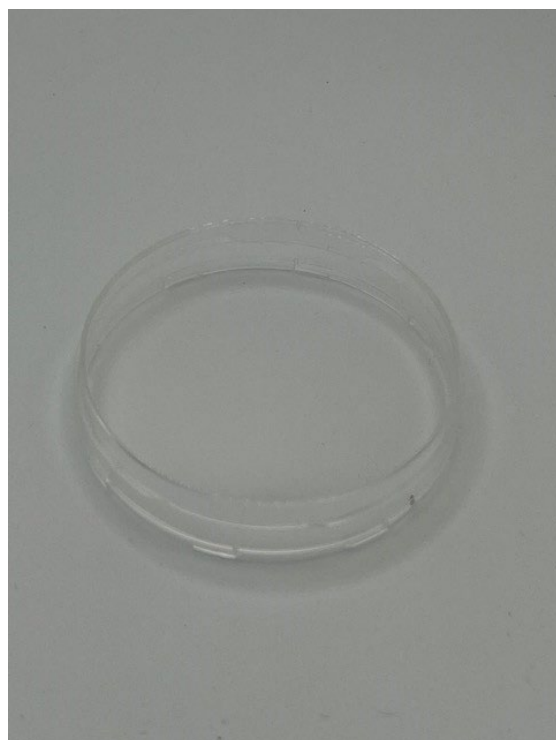



EM Vision lid

EM Vision lid

RMB strives to provide the highest quality products and services to our customers. If you have any questions regarding this informational document, please contact Customer Service at support@rapidmicrobio.com. If interested in placing an order with Rapid Micro Biosystems®, please reach out to orders@rapidmicrobio.com.